

Tomme Tippee®



miomee®
SUPERIOR DESIGN . NATURAL CHOICE

EASY WEIGH BABY SCALES

For further information visit us online
www.tommeetippee.com.au/miomee
www.tommeetippee.co.nz/miomee

Contact us on our Consumer Care Line:

AUS: 1800 890 011

NZ: 0800 877 876



INTRODUCTION

Dear Customer,
Thank you for purchasing the Miomee®
Easy Weigh Baby Scales.

We hope you enjoy using the product.
To ensure your complete satisfaction,
we recommend you read the
instructions carefully.

The Miomee® Easy Weigh Baby Scales
have been designed based on high
quality standards.

The products unique features include;

- **A measurement freeze function**
Allows the weight to be displayed,
even if the baby keeps moving.
- **A highly ergonomic design**
Guarantees comfort and stability.
- **Memory recall and calculation function**
Memorises a measurement and calculates
the difference between the current and
previous weights.

DETAILED PARTS GUIDE



SPECIFICATIONS

- 1 Memorisation of weight (M+)
- 2 Recall of memorised weight (MRC)
- 3 Weight freeze function (HOLD)
- 4 Large 78.5 x 37 mm LCD screen
- 5 On/Off switch (ON/OFF)
- 6 Automatic tare (TARE)

FEATURES

- Maximum weight: 15kg
- Increments: 10g (0-6kgs), 15g (7-15kgs)
- Measurement freeze function
- Weight memory function
- Indicates the difference between the current and previous weights
- Automatic tare
- Power supply: 2 x AAA batteries (included)
- Low battery indicator

OPERATION

1. Putting in the battery

The appliance uses 2 x AAA batteries (included). Insert the 2 x AAA batteries into the battery compartment underneath the baby scale. Respect the polarities +/-.

2. Use

Press the ON/OFF button. All the elements on the screen are displayed for a few seconds, followed by "0.00kg" and the symbol "⌚". When the symbol "⌚" disappears and only "0.00kg" appears on the screen, place the baby in the cradle. If your baby keeps still, his/her weight will automatically stabilise after a few seconds and then appear on the screen. If your baby is moving, the symbol "⌚" appears on the left-hand side of the screen and a non-stabilised weight scrolls on the screen.

Press the HOLD button to freeze the baby's weight. The "⌚" symbol is replaced by "HOLD". The screen displays the baby's stabilised weight. When it is not being used, the baby scale automatically switches off after approximately 20 seconds.

a- Automatic tare function

The automatic tare function allows you to weigh your baby comfortably lying on a blanket without including the weight of the blanket.

- Switch on the scale by pressing the ON/OFF button.
- Put the blanket on the cradle.
- Press the TARE button, the screen displays "---" then "0.00kg".
- Place your baby comfortably on the blanket.

Only the baby's weight is displayed. (When you remove both the baby and the blanket, a negative weight may appear. This is the corresponding weight of the blanket in negative.)

b- Weight memorisation

When the weight of the baby appears on the screen, you can store it in the memory by pressing the "M+" button, and it will continue to be displayed for a few seconds. The "M+" symbol appears.

c- Recalling the recorded weight

Press the "ON/OFF" button to switch on the baby scale. By pressing the "MRC" button, the previous memorised weight will appear on the screen with the sign "Φ" to the right. The memorised weight will continue to be displayed for a few seconds. The recorded weight will remain in the memory until a new weight is recorded and whilst the battery continues working.

d- Difference between two weight measurements

The baby scale automatically displays the difference between the current weight and the previously memorised weight. The difference in weight is displayed using the ▲ or ▼ signs, depending on whether the weight is higher or lower.

Problems	Reason	Solution
Err2	The baby scale is overloaded (weight over 15 kg).	Remove your child immediately or you may damage the baby scale.
ErrE	Baby scale initialisation problem.	Switch the baby scale on and off with no weight on the scale.
bAtt	The battery is worn out.	Replace the battery immediately.

Directions For Use

- Clean your baby scale with a damp cloth without letting water get into the appliance.

! WARNING

TO AVOID INJURY AND DANGER TO BABIES AND CHILDREN, DISPOSE OF ALL PACKAGING IMMEDIATELY

Please read and retain this important information for future reference.

- Remove the two 1.5V batteries if the baby scale is not going to be used for lengthy periods.




Used batteries should be placed in the appropriate collection bin so they can be recycled. Do not mix different types of batteries, or use new and old batteries at the same time. At the end of its serviceable life, this product must be disposed of at a recycling centre for electrical and electronic equipment.

GROWTH RECORD TABLE

NAME	DATE	AGE	TIME	WEIGHT	WEIGHT CHANGE

WORLD HEALTH ORGANISATION (WHO) WEIGHT-FOR-AGE GRAPHS FOR BOYS


These graphs show averages of baby weights from birth to 3 years of age.

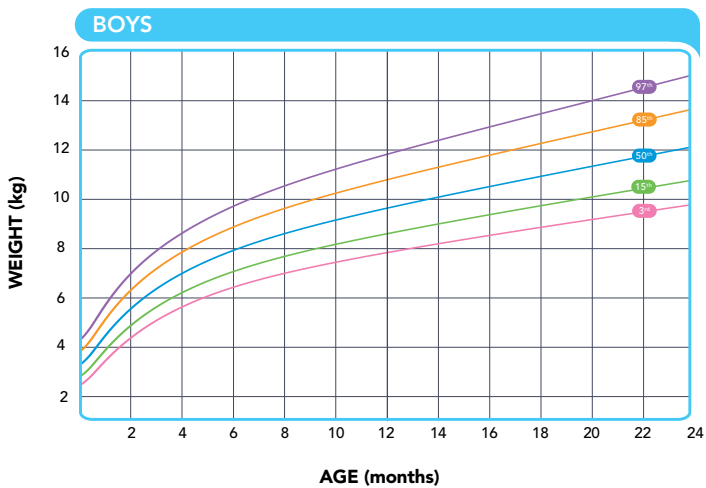

WHO Multicentre Growth Reference Study Group. WHO Child Growth Standards: Length/height-for-age, weight-for-age, weight-for-length, weight-for-height and body mass index-for-age: Methods and development. Geneva: World Health Organization, 2006 (312 pages).

WORLD HEALTH ORGANISATION (WHO) WEIGHT-FOR-AGE GRAPHS

These graphs show averages of baby weights from birth to 3 years of age.

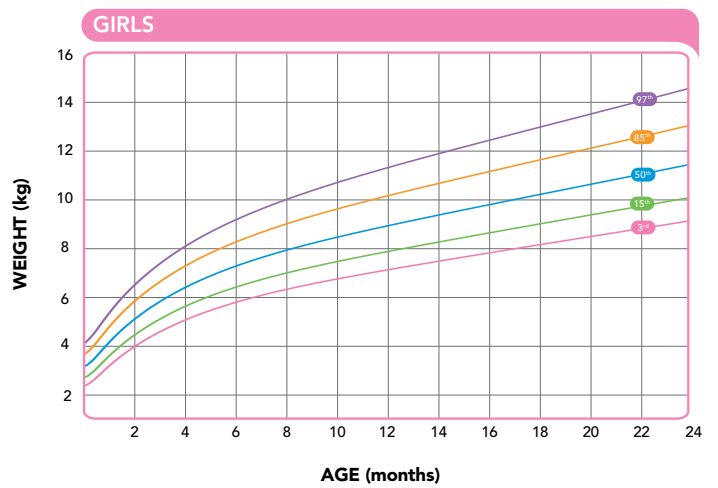

WHO Multicentre Growth Reference Study Group. WHO Child Growth Standards: Length/height-for-age, weight-for-age, weight-for-length, weight-for-height and body mass index-for-age: Methods and development. Geneva: World Health Organization, 2006 (312 pages).

WARRANTY

WARRANTY INFORMATION (AUSTRALIA)

All Miomee® branded products supplied by Jackel Pty Limited ABN 94 840 054 336 to a consumer are guaranteed to be of acceptable quality, subject to the provisions of the Australian Consumer Law. For further information about Australian Consumer Law and consumer guarantees visit www.consumerlaw.gov.au

Any term or condition contained in any warranty given by Jackel Pty Limited in respect of any Miomee® branded product is not intended to reduce or limit any right of a consumer in respect of that product, pursuant to the Trade Practices Act 1974 (Cth), the Competition and Consumer Act 2010 (Cth), or the corresponding legislation of each State and Territory of the Commonwealth of Australia.

Please ensure you keep your original proof of purchase (receipt) and contact our Consumer Care Line on 1800 890 011 or visit www.tommeetippee.com.au/miomee if any problems arise.

WARRANTY INFORMATION (NEW ZEALAND)

There is no requirement to register your purchase of any Miomee® branded product in order to be covered for manufacturer warranty claims. Any term or condition contained in any warranty given by Jackel (NZ) Limited in respect of any Miomee® branded product is not intended to reduce or limit any right of a consumer in respect of that product, pursuant to the Consumer Guarantees Act 1993 or the Fair Trading Act 1986.

Please ensure you keep your original proof of purchase (receipt) and contact our Consumer Care Line on 0800 877 876 or visit www.tommeetippee.co.nz/miomee



Tommee Tippee®



miomee®

SUPERIOR DESIGN . NATURAL CHOICE

are registered Trade Marks of Jackel Trade Marks Pty Limited.

ACN 128 984 690

Copyright© Jackel Trade Marks Pty Limited 2012.

©2012 Imported by Jackel ANZ Pty Limited.

ABN 94 840 054 336

9A, 9 South Street
PO Box 67, Rydalmere
NSW 2116
Australia
Consumer Care Line:
1800 890 011



PO Box 64412
Botany
Manukau City 2163
New Zealand
Consumer Care Line:
0800 877 876

For further information visit us online

www.tommeetippee.com.au/miomee

www.tommeetippee.co.nz/miomee

Product No. 454204

Made in China